

ELATA[®]

E300 Clinical Foam Mattress

WHY WE EXIST

We believe in zero preventable pressure injuries.

It's ambitious but it's necessary. Pressure injuries cause complications. They also add workload and cost. We act before that pain starts.

WHAT THIS MEANS

- **Care** that preserves dignity.
- **Compassion** that simplifies decisions.
- **Commitment** to better outcomes, every day.

HOW WE DELIVER

Advanced systems that adapt to every patient and pathway. Designed to simplify prescription, enhance comfort, and reduce clinical risk.

We deliver personalised, intuitive pressure care, protecting life with dignity, comfort, and trust.



Scan for technical resources

ELATA[®]



ENGINEERED. TESTED. TRUSTED.

WHERE OTHERS SEE ONE FRAME, WE SEE THE WHOLE PICTURE.

We partnered with **Flinders University** to build a **next generation test system**. Traditional mapping captures one static snapshot every 30 minutes. Our system records **thousands of pressure points** every second. This delivers a **true picture** of how weight shifts and movement affects support.

That depth of data underpins Elata's range of mattress models; E300, E400, E500, E600, and E700. So you can trust the results.

AUSTRALIAN DESIGNED & ASSEMBLED. LOCALLY TESTED.

Elata[®] pumps are assembled and rigorously tested here in Australia. That local focus keeps us agile. We listen, adjust, and roll out updates quickly, strengthening performance and delivering the fast, intelligent care clinicians expect.



EVERY PRESSURE INJURY RISK LEVEL. ONE MATCHING SOLUTION.

SIMPLE CHOICES. CLEAR PATHWAYS.

Elata® Range. Quick Category Guide

CATEGORY	MODEL	PURPOSE	PUMP	PG
Home Care	E100, E200	Simple prevention for home and aged care. Default choice in home care pathways.	N/A	9
Clinical Foam	E300	Non-powered. Prevents and treats Stage 1-2 pressure injuries. 9-zone foam with side wall support and heel offloading. Low maintenance.		11
Clinical Hybrid	E400	Use alone for Stage 2 pressure injuries, or step up with a pump to support treatment up to Stage 4. Ideal for home and palliative care. Converts as needs change.	EP10AU	15
Clinical Alternating Hybrid	E500	Alternating air/foam system. Combines comfort with therapy. Treats up to Stage 4 pressure injuries.	EP20AU	19
Clinical Alternating Air	E600	Full alternating air (1-in-2 configuration). Treats up to Stage 4 pressure injuries. Ideal for surgical wards, acute wards, or clients with frailty/reduced mobility.	EP10AU	23
	E700	Advanced alternating (configurable 1-in-4). Auto-setup. Treats up to Stage 4 pressure injuries. For ICU, stable spinal, complex rehab, immobile clients.	EP20AU	27

The EP20AU Pump is fully backwards compatible. It can be used with all Elata pump-ready mattresses (E400 to E700).



Scan to learn more about the range

ELATA[®]

E300

ESSENTIAL COMFORT. INTELLIGENT PREVENTION. WHEN MOST CASES ARE STAGE 1-2 PRESSURE INJURIES. START STRONG.

A high quality specialist **foam mattress** for prevention and treatment. A smart default for home and aged care.

9-Zone Core Foam

Targeted support, total comfort. Precision-engineered into nine zones with varied densities, zoned core foam delivers tailored pressure redistribution across the body's key load points. Each zone adapts to the body to maintain comfort and consistent pressure redistribution.

Accessible Handles

To ensure easy control / transport.

Side Wall Support

Stable transfers; sturdier edges.

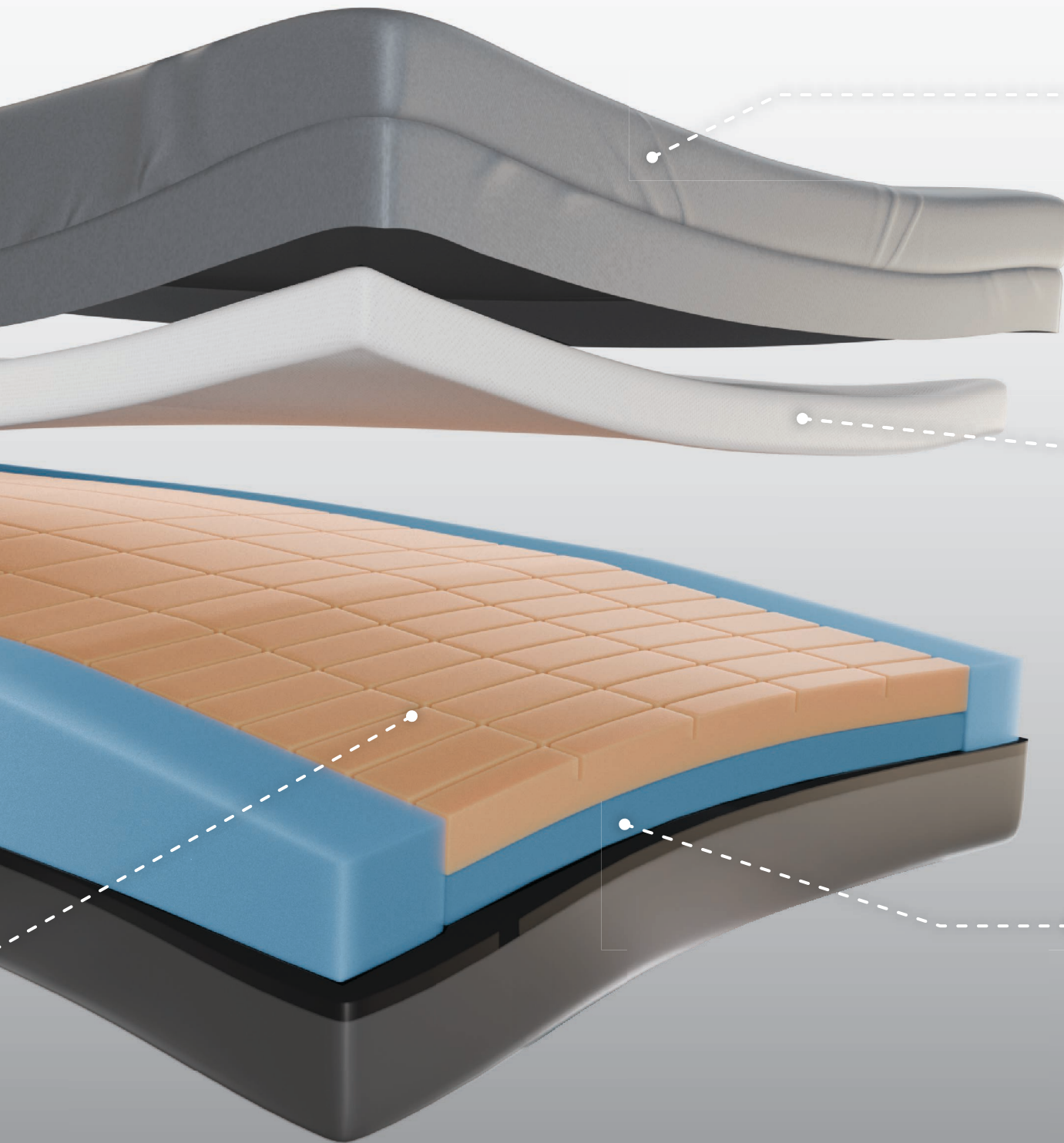
Heel Slope

Gentle decline at the foot end reduces heel pressure and helps prevent injuries.



Cannot be used with a pump
160 kg SWL





Pure Guard Cover

Easy to clean, with a 4-way stretch, vapour-permeable surface that moves with the user to reduce friction and support immersion and envelopment. It allows moisture to escape while blocking fluids, viruses, and bacteria for safer care.

Hydro Lock Barrier

A welded seam and liquid-resistant zip flap prevent fluid penetration, protecting mattress integrity.

Rebound Core Foam

High rebound density and microclimate control. It promotes airflow and delivers lasting pressure redistribution. Maintains skin integrity. Increases ease of mobility.

Anti-Shear

The top layer moves with the client to protect skin and underlying tissues, reducing both visible and hidden shear damage.

Heel Offloading

Foldable foam suspends the heels to relieve 100% of pressure and prevent injuries.

Stability Core Layer

Supportive base layer enhancing stability, comfort, and pressure redistribution.



FOAM MATTRESS E300



Why clinicians choose it:

- 9-Zone Core Foam
- Rebound Core Foam
- Anti-Shear The top layer.
- Heel Slope
- Side Wall Support for stability and edge awareness.

Why it matters:

Better comfort. Fewer changes. Low maintenance. A clear first line of defence.

Most pressure injuries are Stage 1–2 Pressure injuries. E300 helps you stay ahead.

Model E300S	Size: 88 cm * 200 cm * 17 cm SWL: 160 kg
Model E300KS	Size: 107 cm * 200 cm * 17 cm SWL: 160 kg
Warranty	5 years^
Compliance	ARTG: 520915 ISO: 14971, 10993



Scan to learn more about the E300

^Refer to Instructions For Use (IFU)

ELATA

Temperature: 22°C - 71.6°F

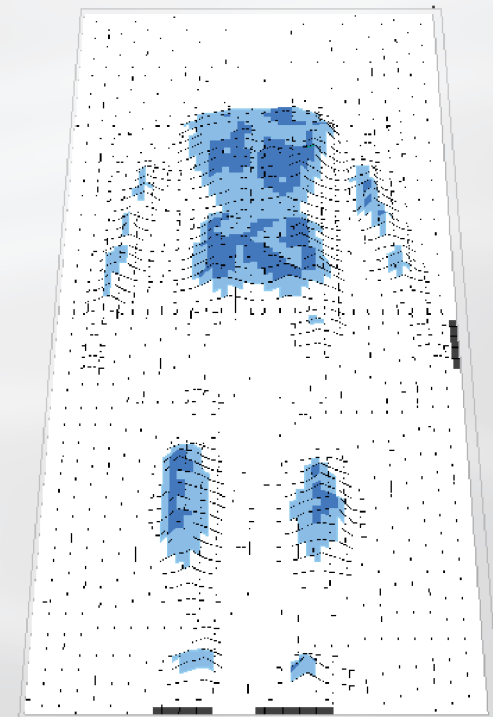
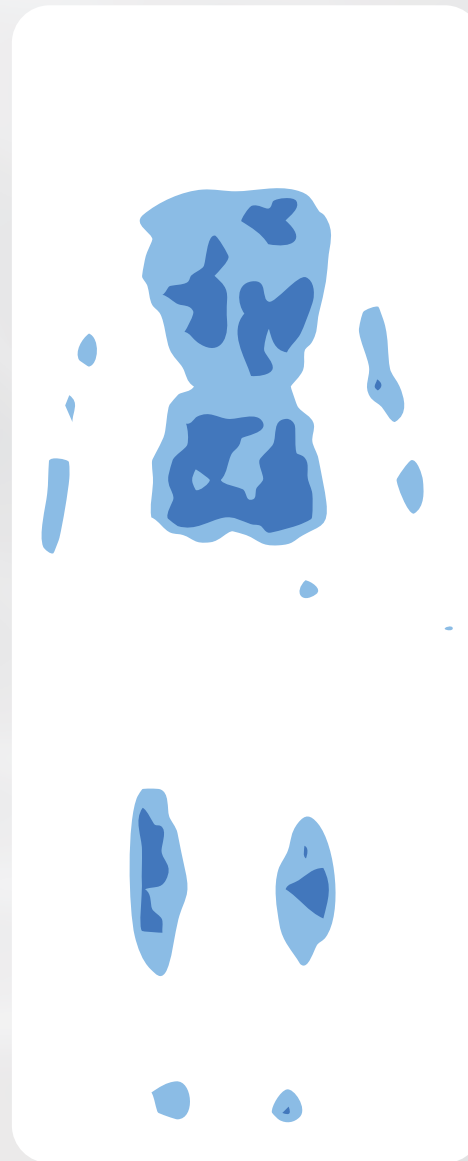
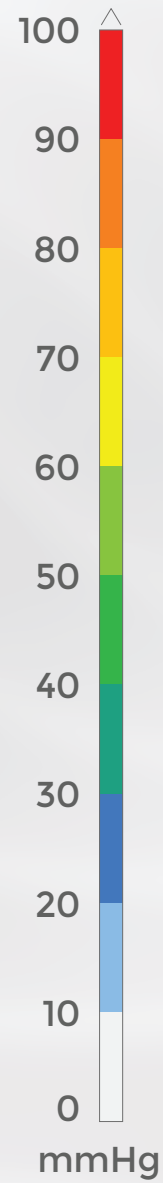
Weight: 75 kg

Duration: 60 mins






Proportional Weighted
95th Percentile Test
Dummy

NON-POWERED

- Mattress: E300
- Date: 2025 / 09 / 25
- MAX: 31.86 mmHg
- Average: 3.20 mmHg
- Frame: 600 / 3600



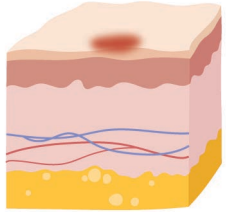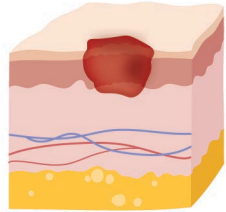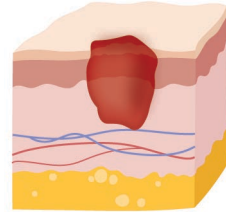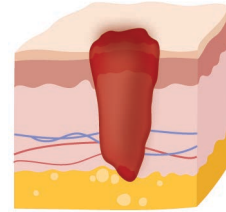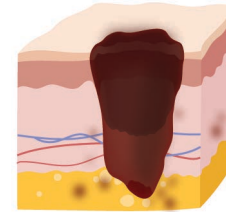PRESSURE INJURY STAGING (NPIAP/EPUAP/PPPIA)

STAGE 1	STAGE 2	STAGE 3	STAGE 4	UNSTAGEABLE
				
<p>Intact skin. Non-blanchable erythema.</p>	<p>Partial-thickness skin loss. Open wound or intact/ruptured serum-filled blister.</p>	<p>Full-thickness skin loss. Fat is visible. Slough may be visible. No exposed bone, muscle or connective tissue.</p>	<p>Full-thickness tissue loss with exposed or palpable bone, muscle or connective tissue. Slough may be visible.</p>	<p>Full-thickness tissue loss, but extent of tissue loss is hidden by slough/eschar.</p>



Scan to book a trial

support@enablelifecare.com.au • 1300 370 370



MATTRESS SELECTOR

THE RIGHT SYSTEM. FIRST TIME.

Visit app.enablelifecare.com.au/assessments, open the Mattress Selector for a guided match based on risk, setting, and clinical goals.



Not sure which system to choose? Download the Enable Lifecare App to quickly match the right Elata® system for every client.

STANDARD IS JUST THE BEGINNING.

// We are passionate about improving pressure care solutions, delivering advanced foam, hybrid, and alternating air mattresses, along with efficient pumps that set new standards in prevention. Through rigorous research and testing, we aim to exceed industry standards. //
We question every detail and ensure each feature serves a purpose.

Ruben Joyce

Ruben Joyce, CEO and Founder

COMPLETE PEACE OF MIND. LONG COVERAGE. REAL PEOPLE. FAST SUPPORT.

- **10 year conditional pump warranty.** Longest in market.
- **In-house service** with fast diagnostics via USB-C log uploads.
- **Education and tools** that simplify trials and prescriptions.

A pressure care mattress is one intervention in a holistic plan. Always assess the whole patient/client.

Scan to see warranty terms



ELATA[®]



support@enablelifecare.com.au

1300 370 370

enablelifecare.com.au